

Dispensette® S

Bottle-top Dispenser



Minimum operating forces,
maximum ease of volume setting.

Proven durability for tough
operating and media conditions.

Number one in
bottle-top dispensing.



Dispensette® S


Innovative ideas with trusted technology – the new bottle-top dispenser Dispensette® S.

- New discharge tube with or without recirculation valve
- New valve system, no sealing rings necessary
- Faster priming due to improved flow technology
- Less force needed during dispensing especially for instruments with large volumes
- Volume selection with interior scalloped track for analog devices, enhances setting reproducibility
- New 1 ml size digital and analog

A Closer Look...

The bottle-top dispenser Dispensette® S has all the features that make dispensing safer and convenient.

Innovative ideas – trusted technology.

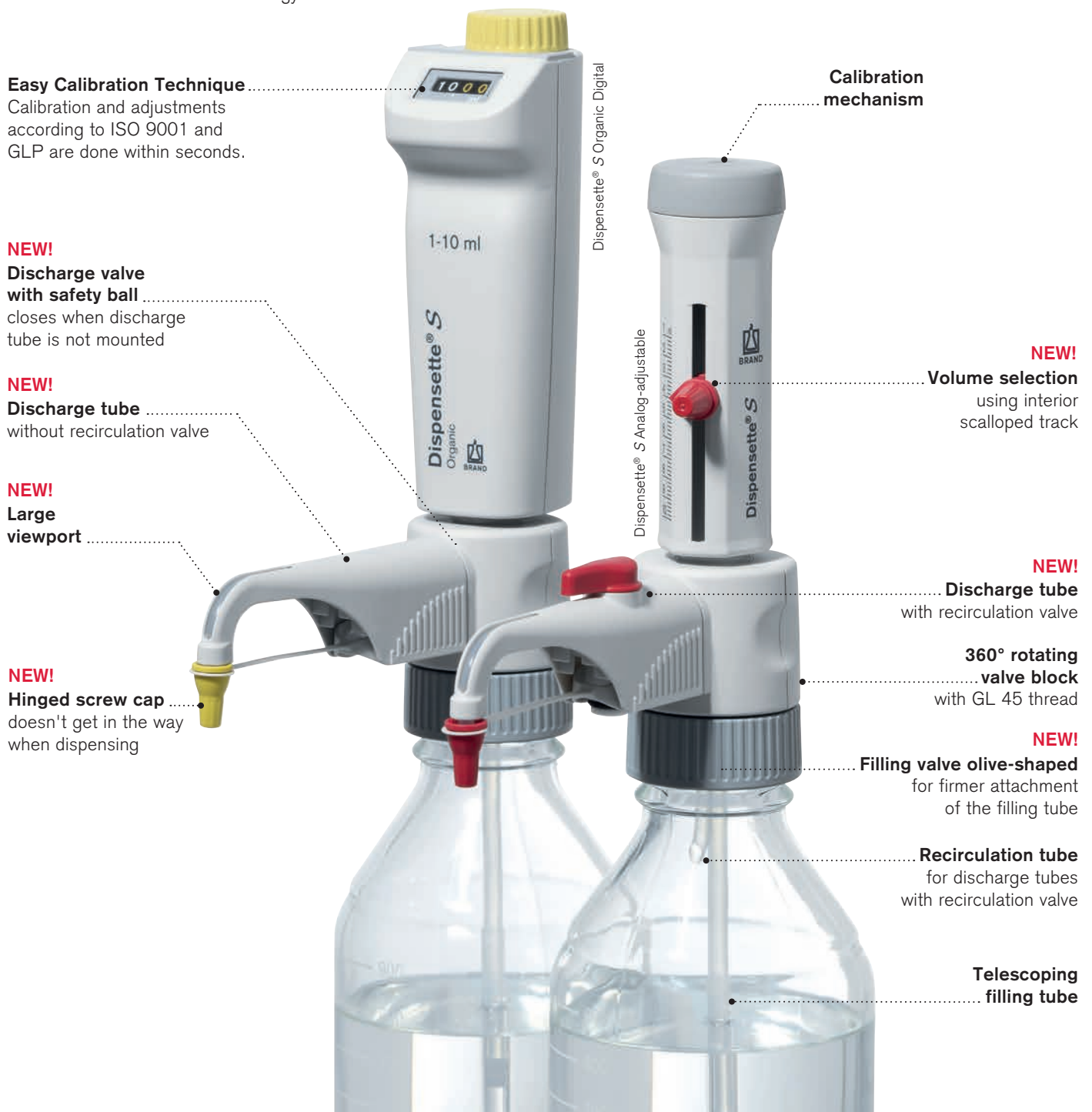

Easy Calibration Technique

Calibration and adjustments according to ISO 9001 and GLP are done within seconds.

NEW!

Discharge valve with safety ball

closes when discharge tube is not mounted

NEW!

Discharge tube

without recirculation valve

NEW!

Large

viewport

NEW!

Hinged screw cap

doesn't get in the way when dispensing

Dispensette® S Organic Digital

Calibration mechanism

Dispensette® S Analog-adjustable

NEW!

Volume selection

using interior scalloped track

NEW!

Discharge tube

with recirculation valve

360° rotating

valve block

with GL 45 thread

NEW!

Filling valve olive-shaped

for firmer attachment of the filling tube

Recirculation tube

for discharge tubes with recirculation valve


Telescoping

filling tube



The right choice

for a wide variety of applications

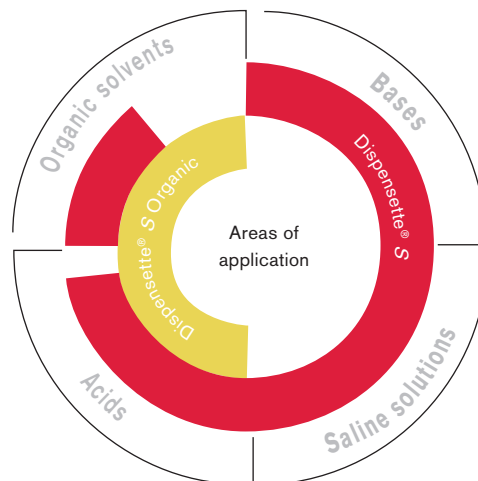

● Dispensette® S

Dispensette® S supports a very wide range of applications for the dispensing of aggressive reagents – directly from the supply bottle:

such as concentrated bases and acids like H_3PO_4 , H_2SO_4 (with certain exceptions such as HCl, HNO_3 , HF, etc.), saline solutions, and a variety of organic solvents.

● Dispensette® S Organic

Dispensette® S Organic is ideal for dispensing organic solvents: such as chlorinated and fluorinated hydrocarbons like trichloroethane and dichloromethane, or acids like concentrated HCl and HNO_3 (except for HF), as well as for trifluoroacetic acid (TFA), tetrahydrofuran (THF), and peroxides.



For dispensing hydrofluoric acid (HF), we recommend the use of the **Dispensette® S Trace Analysis** bottle-top dispenser with platinum-iridium valve spring!

Please find further product information at www.brand.de

Materials in contact with media

● Borosilicate glass, Al_2O_3 -ceramic, platinum-iridium, ETFE, FEP, PFA, PTFE and PP

● Borosilicate glass, Al_2O_3 -ceramic, tantalum, ETFE, FEP, PFA, PTFE and PP

Operating limits

● Vapor pressure max. 600 mbar
viscosity max. 500 mm²/s
temperature max. 40 °C
density max. 2.2 g/cm³

● Vapor pressure max. 600 mbar
viscosity max. 500 mm²/s
temperature max. 40 °C
density max. 2.2 g/cm³



Easy Handling

Easier dispensing and faster priming

In the "floating piston" design, the pistons and cylinders are fitted individually in such a way that the gap between them is just a few micrometers wide. This gap is filled by the fluid being dispensed, forming a zero-wear sealing system with outstanding sliding properties. Optimized flow channels make dispensing and priming even easier – especially in instruments with large volumes.



NEW!
Simple to mount

The new discharge tube is easy to fasten and is available with or without a recirculation valve.



NEW!
Positive volume setting

In analog dispensers, volume setting is quick, secure and repeatable due to the interior scalloped track.



NEW!
Out of the way

So it doesn't impede you while dispensing, the screw cap pivots entirely away from the working area.

Trusted technology



- Autoclavable at 121 °C
- Easy to calibrate and adjust in order to comply with ISO 9001 and GLP guidelines. A positive indicator automatically indicates adjustment from factory settings.
- Easy to dismantle for cleaning
- Replaceable filling and discharge valve with safety ball
- The valve block can be rotated 360° so that the bottle label always faces the user for safety
- Telescoping filling tube adjusts to different bottle sizes
- The 45 mm standard thread plus the included adapters fit common lab bottles
- An extensive line of accessories makes possible special dispensing tasks like sterile applications or dispensing from large containers
- DE-M marking*

* legally replaces  since January 1, 2015

NEW!

Designed without seals

All valves work without any additional sealing rings. That makes cleaning and preparation for autoclaving easier.

Fast calibration

With Easy Calibration technique, you can calibrate quickly and easily in the laboratory in just a few steps. BRAND also offers a factory calibration service.

Serial dispensing

The flexible discharge tube facilitates serial dispensing. It permits fast and precise dispensing even into narrow test tubes.



Dispensing sterile fluids

Dispensette® S Organic and Dispensette® S are completely autoclavable at 121 °C. A connectable microfilter filters the air entering the bottle.



Dispensing sensitive reagents

The drying tube protects sensitive reagents against humidity or CO₂.




Dispenser Selection Chart

Reagent	Dispensette® S		Reagent	Dispensette® S		Reagent	Dispensette® S	
	Organic	Organic		Organic	Organic		Organic	Organic
Acetaldehyde	+	+	Cyclohexane		+	Methylene chloride		+
Acetic acid (glacial), 100%	+	+	Cyclohexanone	+	+	Mineral oil (Engine oil)	+	+
Acetic acid, ≤ 96%	+	+	Cyclopentane		+	Monochloroacetic acid	+	+
Acetic anhydride		+	Decane	+	+	Nitric acid, ≤ 30%	+	+
Acetone	+	+	1-Decanol	+	+	Nitric acid, 30-70% */ **		+
Acetonitrile	+	+	Dibenzyl ether	+	+	Nitrobenzene	+	+
Acetophenone		+	Dichloroacetic acid		+	Oleic acid	+	+
Acetyl chloride		+	Dichlorobenzene	+	+	Oxalic acid	+	+
Acetylacetone	+	+	Dichloroethane		+	n-Pentane		+
Acrylic acid	+	+	Dichloroethylene		+	Peracetic acid		+
Acrylonitrile	+	+	Dichloromethane		+	Perchloric acid	+	+
Adipic acid	+	+	Diesel oil (Heating oil), bp 250-350 °C		+	Perchloroethylene		+
Allyl alcohol	+	+	Diethanolamine	+	+	Petroleum, bp 180-220 °C		+
Aluminium chloride	+	+	Diethyl ether		+	Petroleum ether, bp 40-70 °C		+
Amino acids	+	+	Diethylamine	+	+	Phenol	+	+
Ammonia, ≤ 20%	+	+	Diethylamine	+	+	Phenylethanol	+	+
Ammonia, 20-30%		+	1,2 Diethylbenzene	+	+	Phenylhydrazine	+	+
Ammonium chloride	+	+	Diethylene glycol	+	+	Phosphoric acid, ≤ 85%	+	+
Ammonium fluoride	+	+	Dimethyl sulfoxide (DMSO)	+	+	Phosphoric acid, 85% + Sulfuric acid, 98%, 1:1	+	+
Ammonium sulfate	+	+	Dimethylaniline	+	+	Piperidine	+	+
n-Amyl acetate	+	+	Dimethylformamide (DMF)	+	+	Potassium chloride	+	+
Amyl alcohol (Pentanol)	+	+	1,4 Dioxane	+	+	Potassium dichromate	+	+
Amyl chloride (Chloropentane)		+	Diphenyl ether	+	+	Potassium hydroxide	+	+
Aniline	+	+	Essential oil		+	Potassium permanganate	+	+
Barium chloride	+	+	Ethanol	+	+	Propionic acid	+	+
Benzaldehyde	+	+	Ethanolamine	+	+	Propylene glycol (Propanediol)	+	+
Benzene (Benzol)	+	+	Ethyl acetate	+	+	Pyridine	+	+
Benzine (Petroleum benzin), bp 70-180 °C		+	Ethylbenzene		+	Pyruvic acid	+	+
Benzoyl chloride	+	+	Ethylene chloride		+	Salicylaldehyde	+	+
Benzyl alcohol	+	+	Fluoroacetic acid		+	Scintillation fluid	+	+
Benzylamine	+	+	Formaldehyde, ≤ 40%	+	+	Silver acetate	+	+
Benzylchloride	+	+	Formamide	+	+	Silver nitrate	+	+
Boric acid, ≤ 10%	+	+	Formic acid, ≤ 100%		+	Sodium acetate	+	+
Bromobenzene	+	+	Glycerol	+	+	Sodium chloride	+	+
Bromonaphthalene	+	+	Glycol (Ethylene glycol)	+	+	Sodium dichromate	+	+
Butanediol	+	+	Glycolic acid, ≤ 50%	+	+	Sodium fluoride	+	+
1-Butanol	+	+	Heating oil (Diesel oil), bp 250-350 °C		+	Sodium hydroxide, ≤ 30%	+	+
n-Butyl acetate	+	+	Heptane		+	Sodium hypochlorite	+	+
Butyl methyl ether	+	+	Hexane	+	+	Sulfuric acid, ≤ 98%	+	+
Butylamine	+	+	Hexanoic acid	+	+	Tartaric acid	+	+
Butyric acid	+	+	Hexanol	+	+	Tetrachloroethylene		+
Calcium carbonate	+	+	Hydriodic acid, ≤ 57% **	+	+	Tetrahydrofuran (THF) */ **		+
Calcium chloride	+	+	Hydrobromic acid		+	Tetramethylammonium hydroxide	+	+
Calcium hydroxide	+	+	Hydrochloric acid, ≤ 20%	+	+	Toluene		+
Calcium hypochlorite	+	+	Hydrochloric acid, 20-37% **		+	Trichloroacetic acid		+
Carbon tetrachloride		+	Hydrogen peroxide, ≤ 35%		+	Trichlorobenzene		+
Chloro naphthalene	+	+	Isoamyl alcohol	+	+	Trichloroethane		+
Chloroacetaldehyde, ≤ 45%	+	+	Isobutanol	+	+	Trichloroethylene		+
Chloroacetic acid	+	+	Isooctane		+	Trichlorotrifluoro ethane		+
Chloroacetone	+	+	Isopropanol (2-Propanol)	+	+	Triethanolamine	+	+
Chlorobenzene	+	+	Isopropyl ether	+	+	Triethylene glycol	+	+
Chlorobutane	+	+	Lactic acid	+	+	Trifluoro ethane		+
Chloroform		+	Methanol	+	+	Trifluoroacetic acid (TFA)		+
Chlorosulfonic acid		+	Methoxybenzene	+	+	Turpentine		+
Chromic acid, ≤ 50%	+	+	Methyl benzoate	+	+	Urea	+	+
Chromosulfuric acid	+	+	Methyl butyl ether	+	+	Xylene		+
Copper sulfate	+	+	Methyl ethyl ketone	+	+	Zinc chloride, ≤ 10%	+	+
Cresol		+	Methyl formate	+	+	Zinc sulfate, ≤ 10%	+	+
Cumene (Isopropyl benzene)	+	+	Methyl propyl ketone	+	+			

The above recommendations reflect testing completed prior to publication. Always follow instructions in the operating manual of the instrument as well as the reagent manufacturer's specifications. In addition to these chemicals, a variety of organic and inorganic saline solutions (e.g., biological buffers), biological detergents and media for cell culture can be dispensed. Should you require information on chemicals not listed, please feel free to contact BRAND. Status as of: 1116/13

* use ETFE/PTFE bottle adapter
** use PTFE seal for valve block

Note!  For dispensing HF, we recommend the use of the Dispensette® S Trace Analysis bottle-top dispenser with platinum-iridium valve spring.



Ordering Information

Items supplied:

Dispensette® S / Dispensette® S Organic bottle-top dispenser, DE-M marking, performance certificate, telescoping filling tube, recirculation tube (optional), mounting tool and adapters of polypropylene:

Nominal volume ml	Adapter for bottle thread	Filling tube length
1, 2, 5, 10	GL 24-25, GL 28/S 28, GL 32-33, GL 38, S 40	125-240 mm
25, 50, 100	GL 32-33, GL 38, S 40	170-330 mm

Dispensette® S

Capacity ml	Subdivision ml	A* ≤ ±		CV* ≤		without recirculation valve Cat. No.	with recirculation valve Cat. No.
		%	µl	%	µl		
■ Dispensette® S, Digital							
0.1 - 1	0.005	0.6	6	0.2	2	4600 310	4600 311
0.2 - 2	0.01	0.5	10	0.1	2	4600 320	4600 321
0.5 - 5	0.02	0.5	25	0.1	5	4600 330	4600 331
1 - 10	0.05	0.5	50	0.1	10	4600 340	4600 341
2.5 - 25	0.1	0.5	125	0.1	25	4600 350	4600 351
5 - 50	0.2	0.5	250	0.1	50	4600 360	4600 361
■ Dispensette® S, Analog-adjustable							
0.1 - 1	0.02	0.6	6	0.2	2	4600 100	4600 101
0.2 - 2	0.05	0.5	10	0.1	2	4600 120	4600 121
0.5 - 5	0.1	0.5	25	0.1	5	4600 130	4600 131
1 - 10	0.2	0.5	50	0.1	10	4600 140	4600 141
2.5 - 25	0.5	0.5	125	0.1	25	4600 150	4600 151
5 - 50	1.0	0.5	250	0.1	50	4600 160	4600 161
10 - 100	1.0	0.5	500	0.1	100	4600 170	4600 171
■ Dispensette® S, Fixed-volume							
1		0.6	6	0.2	2	4600 210	4600 211
2		0.5	10	0.1	2	4600 220	4600 221
5		0.5	25	0.1	5	4600 230	4600 231
10		0.5	50	0.1	10	4600 240	4600 241
Special fixed volumes: 0.5-100 ml (please state when ordering)						4600 290	4600 291




Dispensette® S Organic

Capacity ml	Subdivision ml	A* ≤ ±		CV* ≤		without recirculation valve Cat. No.	with recirculation valve Cat. No.
		%	µl	%	µl		
■ Dispensette® S Organic, Digital							
0.5 - 5	0.02	0.5	25	0.1	5	4630 330	4630 331
1 - 10	0.05	0.5	50	0.1	10	4630 340	4630 341
2.5 - 25	0.1	0.5	125	0.1	25	4630 350	4630 351
5 - 50	0.2	0.5	250	0.1	50	4630 360	4630 361
■ Dispensette® S Organic, Analog-adjustable							
0.5 - 5	0.1	0.5	25	0.1	5	4630 130	4630 131
1 - 10	0.2	0.5	50	0.1	10	4630 140	4630 141
2.5 - 25	0.5	0.5	125	0.1	25	4630 150	4630 151
5 - 50	1.0	0.5	250	0.1	50	4630 160	4630 161
10 - 100	1.0	0.5	500	0.1	100	4630 170	4630 171
■ Dispensette® S Organic, Fixed-volume							
5		0.5	25	0.1	5	4630 230	4630 231
10		0.5	50	0.1	10	4630 240	4630 241
Special fixed volumes: 2-100 ml (please state when ordering)						4630 290	4630 291



* Calibrated to deliver (TD, Ex). Error limits according to the nominal capacity (= maximum volume) indicated on the instrument, obtained with instrument and distilled water at equilibrium with ambient temperature at 20 °C, and with smooth, steady operation. The error limits are well within the limits of DIN EN ISO 8655-5. DE-M marking. A = Accuracy, CV = Coefficient of variation

Note!  For trace analysis and dispensing HF, we recommend the use of the Dispensette® S Trace Analysis bottle-top dispenser.

Accessories · Spare Parts



Discharge tubes

With and without recirculation valve. Screw cap PP.
Pack of 1.

Description	Nominal volume ml	Shape	Length mm	without recirculation valve Cat. No.	with recirculation valve Cat. No.
Dispensette® S	1, 2, 5, 10	fine tip	108	7080 02	7081 02
	5, 10	standard	108	7080 05	7081 04
	25, 50, 100	fine tip	135	7080 06	7081 06
	25, 50, 100	standard	135	7080 08	7081 09
Dispensette® S Organic	5, 10	standard	108	7080 14	7081 14
	25, 50, 100	fine tip	135	7080 16	7081 16
	25, 50, 100	standard	135	7080 19	7081 19

Telescoping filling tubes

For Dispensette® S and Dispensette® S Organic.
FEP. Adjusts to various bottle heights.
Pack of 1.



Nominal volume ml	Outer Ø mm	Length mm	Cat. No.
1, 2, 5, 10	6	70-140	7082 10
		125-240	7082 12
		195-350	7082 14
		250-480	7082 16
25, 50, 100	7.6	170-330	7082 18
		250-480	7082 20

Flexible discharge tube with recirculation valve*

For Dispensette® S and Dispensette® S Organic.
PTFE, coiled, length approx. 800 mm, with safety handle.
Pack of 1.



Nominal volume ml	Discharge tube Outer Ø mm	Inner Ø mm	Cat. No.
2, 5, 10	3	2	7081 32
25, 50, 100	4.5	3	7081 34

* not suitable for HF

Bottle stand

PP. Full plastic construction.
Support rod 325 mm,
base plate 220 x 160 mm,
weight 1130 g.
Pack of 1.



Cat. No. 7042 75

Sealing ring for valve block


PTFE. For highly volatile media.
Pack of 1.



Cat. No. 7044 86

Drying tube incl. PTFE-sealing ring

Without drying agent.
Pack of 1.

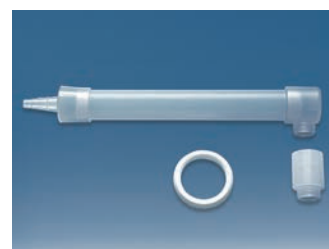

Cat. No. 7079 30

Additional accessories can be found at www.brand.de

BRAND®, Dispensette®, BRAND. For lab. For life.® and the word-image trademarks here are brands of BRAND GMBH + CO KG, Germany.

Our technical literature is intended to inform and advise our customers. However, the validity of general empirical values, and of results obtained under test conditions, for specific applications depends on many factors beyond our control. Please appreciate, therefore, that no claims can be derived from our advice. The user is responsible for checking the appropriateness of the product for any particular application.

Subject to technical modification without notice. Errors excepted.

BRAND GMBH + CO KG · P.O. Box 11 55 · 97861 Wertheim · Germany
Tel.: +49 9342 808-0 · Fax: +49 9342 808-98000 · E-Mail: info@brand.de · Internet: www.brand.de